
ALGILOCK®

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Impression tray

Porte-empreinte

Avtryckssked

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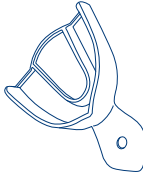EN

Processing instructions according to EN ISO 17664:2021

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Warnings: Caution! The product must be cleaned, disinfected and sterilized according to these processing instructions before its first use as well as after every subsequent use. Disinfection alone is not sufficient. Thorough cleaning and disinfection are decisive prerequisites for effective sterilization. The processing should be started as soon as possible but not later than 2 hours after use. During the processing, the product should not be unnecessarily exposed to moisture or humidity. Please also observe the legal regulations applicable in your country as well as the hygiene regulations of the medical practice and/or the hospital. **Limitations in the processing:** Frequent processing has only minor effects on this product. The end of the product's useful life is primarily determined by wear and damage caused by use. In case of doubt, questionable products should always be sorted out and replaced early. The user alone is responsible for making the decision regarding repeated use of the product. If the product is used too frequently, the manufacturer does not accept any warranty for the function, performance and safety of the product. **Personal protective equipment:** For reasons of occupational safety and to minimize infections/cross infections, suitable personal protective equipment (protective clothes, protective gloves, protective goggles and nose/mouth protection) must be used during the entire processing.

INSTRUCTIONS:

INITIAL TREATMENT AT THE PLACE OF USE

Remove surface contamination using a disposable cloth/paper cloth while still at the place of use. Rinse

the products with water (minimum drinking water quality) no later than 2 hours after use. Do not allow residues or contamination of any kind to dry on the product. When using the product for the first time at the place of use, do not use any agents containing aldehyde or alcohol as they may cause protein fixation.

PREPARATION BEFORE CLEANING

If possible, disassemble the product into its individual components. Visually inspect for damage and wear.

CLEANING & DISINFECTION: BY MACHINE (RECOMMENDED)

Only perfectly maintained, checked, calibrated and released cleaning and disinfection equipment according to EN ISO 15883 with valid test mark (CE mark) may be used. When cleaning and disinfecting several different products/individual components by machine, ensure that the cleaning device is not overfilled and that the products/individual components do not come into contact with each other. We recommend completing the cleaning and disinfection cycle according to the specifications of the device manufacturer and the manufacturer of the cleaning agent or disinfectant. Ensure that the cleaning agent is compatible with the products and disinfectants that might be used. The pH of the cleaning agent and disinfectant should be between 5.5 and 8.5. Do not use organic solvents (e.g. alcohols, ether, ketones, petrol), oxidizing agents (e.g. peroxides), halogens (chlorine, iodine, bromine) or aromatic/halogen hydrocarbons. We recommend the following products: Dr. Weigert, neodisher MA/Dr. Weigert, neodisher Z. The manufacturer's instructions for use of the cleaning agent and disinfectant must strictly be observed. In particular, note the concentrations and exposure times to be used.

Description of the recommended cleaning process:

Pre-cleaning: The pre-cleaning should be completed in a water bath (minimum drinking water quality) without other cleaning additions using a cleaning brush below the water surface to achieve both sufficient cleaning of the products/individual components while avoiding protein fixation, and to protect the environment from contamination with splash water. The duration should be based on the level of contamination of the product/individual component but should not be less than 1 minute. Information on the cleaning and disinfection device: Type: Miele G7881, Device for pre-cleaning: Flexbrush, REF 605 254

Description of the recommended program:

Cleaning program: VarioTD **Cleaning temperature:** 50–56.28 °C **Hold time:** 10:56 min. **Disinfection temperature:** 65–94.41 °C **Hold time:** 09:59 min. **Flushing technology:** Fresh water system/demineralized water **Cleaning agent:** Dr. Weigert, neodisher MA, **Dosing:** 30 mg (3 g/L) **Neutralisator:** Dr. Weigert, neodisher Z, **Dosing:** 1–2 ml/L **Aids:** Insert rack for impression trays
Use only water that is drinking water quality. Rinsing and drying should be completed by machine. The machine cleaning should not exceed a temperature of 60 °C. Ensure that visible contamination has been completely removed. If visible contamination remains after the cleaning and disinfection process, the cleaning & disinfection must be repeated. Before further processing, the product must be free from any residues and dry.

CLEANING: MANUAL

Manual cleaning should be completed in a water bath (minimum drinking water quality) with the cleaning agent specified below using a brush below the water surface to achieve both sufficient cleaning of the products/individual components while avoiding protein fixation, and to protect the environment from contamination with splash water. The duration should be based on the level of contamination of the product/individual component but should not be less than 1 minute. Subsequently, rinse the products/individual components with water (minimum drinking water quality) for at least 20 seconds.

Description of the recommended cleaning process:

Pre-cleaning: Remove surface contamination using a disposable cloth **Cleaning:** Insert the instrument into the strainer bowl. Hang the strainer bowl into the cleaning bath with cleaning solution, remove contamination using cleaning brushes. **Cleaning agent:** ID 213, company Dürr Dental, concentration: 2 % **Rinsing:** Rinse with water (minimum drinking water quality) for 20 sec. **Drying:** Dry at room temperature
Manual cleaning should not exceed a temperature of 45 °C. Ensure that visible contamination has been completely removed. If visible contamination remains after the cleaning process, the cleaning must be repeated.

DISINFECTION: MANUAL

For manual disinfection, please use only approved disinfectants with tested effectiveness (CE mark, VAH-/DGHM-listed). Put the products/individual components into corresponding disinfection baths according to the specifications of the disinfectant manufacturer. Ensure that the products are sufficiently covered and do not contact each other. The instructions for use of the disinfectant manufacturer must be strictly observed. In particular, note the concentrations and exposure times to be used. We recommend the following products: Dürr Dental ID 213 instrument disinfection

Description of the recommended disinfection process:

Carrying out proper cleaning preparation and manual cleaning

Disinfection: Insert the instrument into the strainer bowl. Hang the strainer bowl into the disinfection bath with disinfection solution **Disinfectant:** ID 213, company Dürr Dental, concentration: 2 % **Contact time:** 5 min. **Rinsing:** Rinse with water (minimum drinking water quality) for 20 sec. **Drying:** Dry at room temperature

Ensure that the disinfectant is compatible with the products and cleaning agents that might be used. The pH of the disinfectant should be between 5.5 and 8.5. Do not use organic solvents (e.g. alcohols, ether, ketones, petrol), oxidizing agents (e.g. peroxides), halogens (chlorine, iodine, bromine) or aromatic/halogen hydrocarbons. Please observe thermal restrictions that might apply, taking into account the disinfectant manufacturer's specifications. Subsequently rinse the products/individual components with water (minimum drinking water quality) for at least 20 seconds. Before further processing, the product must be free from any residues and dry.

DRYING

If you process the products manually, the products/individual components may be dried by blowing them off with filtered, oil-free compressed air according to DIN ISO 8573-1 (medical cleanliness class) or at room temperature. Drying in the machine reprocessing method should take place directly through devices with hot air. If machine drying, do not exceed a temperature of 100 °C.

MAINTENANCE, CONTROL AND TESTING

There are no special maintenance requirements for the product. Re-assemble disassembled products/individual components. A visual inspection for contamination, damage, wear, deformation must always take place before and after the individual work steps. Damaged or corroded products must no longer be used. If the product/individual component is not visibly clean, the entire processing must be repeated or the product/individual component must be disposed of properly.

PACKAGING

The product must be packed in suitable and standardized transparent sterile packaging (sterilization bags) and sealed. Observe the instructions of the sterilization bag and sealing machine manufacturers as well as the current, normative requirements. Products/individual parts that have not been sterilized in bags must be used immediately.

STERILIZATION

Use only tested steam vacuum autoclaves. When sterilising several different products/individual components, ensure that the autoclave is not overfilled and that the products/individual components are not in contact with each other. The following sterilization cycles may be completed: Steam sterilization, 134 °C, hold time 5 minutes or steam sterilization, 121 °C, hold time 15 minutes. To dry the products/individual components, the drying cycle of the autoclave should be set. Use the autoclave manufacturer's instructions for use.

STORAGE

To maintain the sterility, store the products in standardized sterilization bags at a dry, clean place until they are used. If the sterile packaging is damaged, reprocess the products again before use.

ADDITIONAL INFORMATION

Label the sterilized products according to the legal and national regulations. The recommended storage period for sterile medical devices is described in standard DIN 58953-8 and depends on external influences and effects during storage, transport and handling. The user must ensure within the scope of their quality management system that specified processing cycles that might apply (see Limitations in the processing) are not exceeded. The manufacturer and the competent authority of your member state must be immediately notified all serious incidents occurring in connection with the product. Recommendations and regulations for the disposal of the products vary by region. Ask your competent disposal company for the current regulations in your region.

MANUFACTURER CONTACT DETAILS

Phone: (+49) 203/99269-0, Fax: (+49) 203/299283, Email: info@hagerwerken.de

The instructions specified above are validated by the medical device manufacturer as suitable for the preparation of a medical device for its reuse. The processor is responsible for ensuring that the reprocessing actually carried out with the equipment used, materials and personnel in the preparation facility achieves the desired result. Verification and/or validation and routine monitoring of the process are required for this purpose.

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HAGER & WERKEN GmbH & Co. KG · Ackerstraße 1
Tecor-Dental AG · Ebnetstrasse 65 8200 Schaffhausen, Schweiz