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
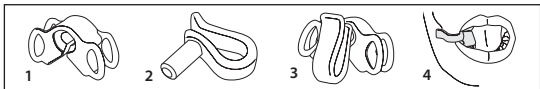
MIRAHOLD®-BLOCK

EN Instructions for Use

- Elastic bite blocks for mouth opening
- Silicone block is positioned between upper and lower jaw and keeps the mouth reliably open by spring force
- Fixation possible with cheek retractor
- Autoclavable

Bite-Block Mirahold-Block for mouth opening with bite block and retractor. The bite block is used as shown in pictures 1 – 4.

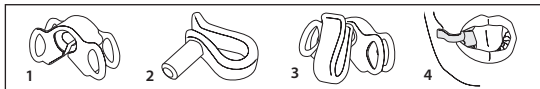Processing instructions according to EN ISO 17664:2021

Product: Mirahold-Block

REF 605 210 **REF** 605 238 **REF** 605 220 **REF** 605 239

Manufacturer: HAGER & WERKEN GmbH & Co. KG

Ackerstraße 1, 47269 Duisburg, Germany


Warnings: Attention! The product must be cleaned, disinfected and sterilized according to these processing instructions before its first use as well as after every further use. Disinfection alone is not enough. Thorough cleaning and disinfection are essential prerequisites for effective sterilization. The processing should be started as early as possible, at latest 2 hours after use. During the processing, the product should not be unnecessarily exposed to wetness or humidity. Please also observe the legal regulations applicable in your country as well as the hygiene regulations of your medical practice and/or the hospital. **Limitations in the processing:** Frequent reprocessing has only minor effects on this product. The end of the product's useful life is mainly determined by wear and damage caused by use. In case of doubt, the products should always be sorted out and replaced at an early stage. The user is solely responsible for making the decision regarding repeated use of the product. If the product is used too frequently, the manufacturer does not accept any warranty for the function, performance and safety of the product. **Personal protective equipment:** For reasons of occupational safety and to minimize infections/cross-infections, suitable personal protective equipment (protective clothes, protective gloves, protective goggles and nose/mouth protection) must be used during the entire processing procedure.

INITIAL TREATMENT AT THE PLACE OF USE

Remove surface contamination using a disposable cloth/paper cloth while still at the place of use. Rinse the products with water (at least drinking water quality) 2 hours after use at the latest. The drying of residues or contamination of any kind on the product must be avoided. For the initial treatment at the place of use, do not use any aldehyde-containing or alcohol-containing agents as they may cause protein fixation.

PREPARATION BEFORE CLEANING

Disassemble the product into its individual components. Visual inspection for damage and wear.

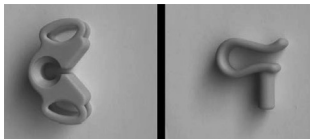CLEANING & DISINFECTION: BY MACHINE (RECOMMENDED)

Only perfectly maintained, checked, calibrated and approved cleaning and disinfection equipment according to EN ISO 15883 with a valid check mark (CE mark) may be used. During machine cleaning and disinfection of multiple different products/individual components, make sure that the cleaning device is not overfilled and that the products/individual components do not come into contact with each other. We recommend completing the cleaning and disinfection cycle according to the device manufacturer's and the cleaning agent and/or disinfectant manufacturer's specifications. Make sure that the cleaning agent is compatible with the products and any disinfectants that are used. The pH of the cleaning agent and disinfectant should lie between 5.5 and 8.5. Do not use organic solvents (e.g. alcohols, ethers, ketones, benzines), oxidizing agents (e.g. peroxides), halogens (chlorine, iodine, bromine) or aromatic/halogen hydrocarbons. Use with the following agents is recommended: - Dr. Weigert, neodisher MA/- Dr. Weigert, neodisher Z. The usage instructions of the device manufacturer and the manufacturer of the cleaning agent and disinfectant must be strictly observed. In particular, concentrations and contact times must be observed.

Description of the recommended cleaning process: Pre-cleaning: The pre-cleaning should be completed in a water bath (at least drinking water quality) without other cleaning additives using a cleaning brush below the water surface to achieve sufficient cleaning of the products/individual components while avoiding protein fixation and also to protect the environment from contamination with splash water. The duration should be based on the level of contamination of the product/individual component; it should not however be less than 1 minute.

Cleaning tools for pre-cleaning: Mirasuc Brush, REF 605 216/Chiru-Cleaner 23, REF 605 215

Information about the cleaning and disinfection device: Type: Miele G7881

Description of the recommended cleaning and disinfection:

Cleaning program: VarioTD

Cleaning temperature: 50-56.6°C

Hold time: 09:55 min

Disinfection temperature: 65 – 94.4°C

Hold time: 09:55 min

Rinsing technique: Fresh water system/demineralized water

Cleaning agent: Dr. Weigert, neodisher MA, dosing: 30 g (3 g/l)

Neutralizer: Dr. Weigert, neodisher Z, dosing: 1 – 2 ml/l.

The water should be of drinking water quality. Rinsing and drying should be completed by machine. Machine cleaning should not exceed a temperature of 60°C. It must be ensured that visual contamination has been completely removed. If there is still visible contamination after the cleaning and disinfection process, the cleaning & disinfection must be repeated. Before further processing, the product must be free from any residues and dry.

CLEANING: MANUAL

Manual cleaning should be completed in a water bath (at least drinking water quality) with the cleaning agent specified below using a brush below the water surface to achieve sufficient cleaning of the products/individual components while avoiding protein fixation and also to protect the environment from contamination with splash water. The duration should be based on the level of contamination of the product/individual component; it should not however be less than 1 minute. Afterwards, rinse the products/individual components with water (at least drinking water quality) for at least 20 seconds.

Description of the recommended cleaning process:

Pre-cleaning: Remove surface contamination using a disposable cloth

Cleaning: Insert the instrument into the strainer bowl. Suspend the strainer bowl in the cleaning bath with cleaning solution, remove contamination by means of cleaning brushes.

Cleaning agent: ID 213, made by Dürr Dental, concentration: 2%

Cleaning tools: Mirasuc Brush, REF 605 216/Chiru-Cleaner 23, REF 605 215

Post-rinsing: Rinse with water (at least drinking water quality) for 20 seconds

Drying: Dry at room temperature

Manual cleaning should not exceed a temperature of 45°C. It must be ensured that visual contamination has been completely removed. If there is still visible contamination after the cleaning process, cleaning must be repeated.

DISINFECTION: MANUAL

For manual disinfection, please only use approved disinfectants with tested effectiveness (CE mark, VAH-/DGHM-listed). Put the products/individual components into suitable disinfection baths according to the disinfectant manufacturer's specifications. Ensure that the products are sufficiently covered and do not come into contact with each other. The usage instructions of the disinfectant manufacturer must be strictly observed. In particular, concentrations and contact times must be observed.

Use of the following agent is recommended: - Dürr Dental ID 213 instrument disinfection

Description of the recommended disinfection process: Proper cleaning preparation and manual cleaning

Disinfection: Insert the instrument into the strainer bowl. Suspend the strainer bowl in the disinfection bath with disinfectant solution

Disinfectant: ID 213, made by Dürr Dental, concentration: 2%

Cleaning tools: Mirasuc Brush, REF 605 216/Chiru-Cleaner 23, REF 605 215

Contact time: 5 min

Post-rinsing: Rinse with water (at least drinking water quality) for 20 seconds

Drying: Dry at room temperature

Make sure that the disinfectant is compatible with the products and any cleaning agents that might be used. The pH of the disinfectant should lie between 5.5 and 8.5. Do not use organic solvents (e.g. alcohols, ethers, ketones, benzines), oxidizing agents (e.g. peroxides), halogens (chlorine, iodine, bromine) or aromatic/halogen hydrocarbons. Please observe any thermal restrictions that might apply while taking the disinfectant manufacturer's specifications into consideration. Afterwards, rinse the products/individual components with water (at least drinking water quality) for at least 20 seconds. Before further processing, the product must be free from any residues and dry.

DRYING

If you process the products manually, the products/individual components may be dried by blowing with filtered, oil-free compressed air according to DIN ISO 8573-1 (medical purity class) or at room temperature. If you process the products with machines, drying should be performed directly by the relevant equipment using hot air. For machine drying, a temperature of 100°C should not be exceeded.

MAINTENANCE, INSPECTION AND TESTING

The product does not have any special maintenance requirements. Reassemble any disassembled products/individual components. You must always carry out a visual inspection for contamination, damage, wear and deformation before and after each work step. Damaged or corroded products must no longer be used. If the product/individual component is not visibly clean, the full processing procedure must be repeated or the product/individual component must be properly disposed of.

PACKAGING

The product must be packed in suitable and standardized transparent sterile packaging (sterilization bags) and sealed. Observe the instructions of the manufacturers of the sterilization bag and sealing machine as well as the current normative requirements. Products/individual components that have not been sterilized in bags must be used immediately.

STERILIZATION

Only tested steam vacuum autoclaves may be used. During the sterilization of multiple different products/individual components, make sure that the autoclave is not overfilled and that the products/individual components do not come into contact with each other. The following sterilization cycles may be completed: Steam sterilization, 134°C, hold time 5 minutes; or steam sterilization, 121°C, hold time 20 minutes. To dry the products/indi-

vidual components, the drying cycle of the autoclave should be configured accordingly. Please observe the autoclave manufacturer's usage instructions.

STORAGE

To maintain sterility, products must be stored in standardized sterilization bags in a dry, clean place until they are used. If the sterile packaging is damaged, the products must be processed again before they are used.

ADDITIONAL INFORMATION

Mark the sterilized products according to the legal and national regulations. The recommended storage duration for sterile medical devices is described in the DIN 58953-8 standard and depends on external influences and effects during storage, transport and handling. The user must ensure within their quality management system that any specified processing cycles that might apply (see Limitations in the processing) are not exceeded. The manufacturer and the competent authority of your member state must be immediately notified of all serious incidents occurring in connection with the product. For disposal of products, recommendations and regulations differ according to region. Ask your disposal company for the current regulations in your region.

CONTACT TO THE MANUFACTURER

Tel.: (+49) 203/99269-0 Fax.: (+49) 203/299283 Email: info@hagerwerken.de

The instructions specified above are validated by the medical device manufacturer as suitable for the preparation of a medical device for its reuse. The processor is responsible for ensuring that the processing performed in practice, including the equipment used, materials and personnel in the preparation facility, achieves the desired result. To this end, verification and/or validation and routine monitoring of the procedure are required.

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